




## UA-787 Quick Response with EasyCuff™ Feature

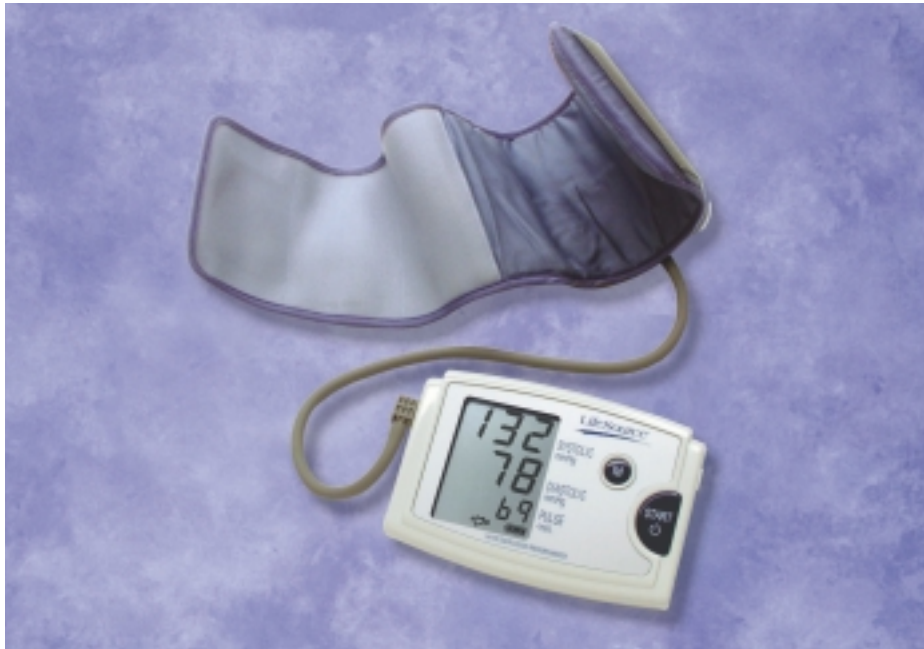

The EasyCuff and Irregular Heartbeat features allow this Blood Pressure Monitor to adapt to more users than ever before.

The Quick Response with *EasyCuff* feature gives fast and accurate blood pressure readings on a large, easy to read display screen while retaining the last 30 blood pressure and pulse rate measurements in memory.

Included with the blood pressure monitor is the *EasyCuff* feature, the latest innovation in cuff design from A&D Medical. With its unique semi-hard shell construction, the cuff readily conforms to the shape of your upper arm for a quick and proper fit.

The Quick Response monitor includes a unique Irregular Heartbeat feature. The UA-787 monitor now provides blood pressure and pulse rate measurements even when an Irregular Heartbeat occurs.

- Simultaneous readout of systolic/diastolic pressure and pulse rate
- Stores last 30 blood pressure and pulse rate measurements in memory
- *EasyCuff* size is two inches larger than the standard medium cuff to fit more arm sizes
- Large, easy to read 20mm digital display
- Large-print instruction guide with illustrations and quick start card in English, French and Spanish
- Four week blood pressure logbook included

For questions about *LifeSource* monitors, call our toll free *LifeSource* HealthLine at 1-888-726-9966 and a specially trained representative will assist you, or visit us at [www.LifeSourceOnline.com](http://www.LifeSourceOnline.com).

- Professional accuracy
- Fast measurement
- Irregular Heartbeat feature
- 30 reading memory
- Includes innovative easy to use cuff
- Extra large digital display
- One button operation
- Quiet operation
- Latex free



UA-787R



UA-787V



# UA-787 Quick Response with EasyCuff™ Feature

To obtain an accurate reading, it is important to use the correct cuff size. A cuff that is too small will result in a higher than expected reading and a cuff that is too large will result in a lower than expected reading. The **EasyCuff** feature has a larger cuff range (7.5" to 14.2" (19-36cm)) than the size of the standard medium cuff (7.5" to 12.2" (19-31cm)) to fit more arm sizes.

UNITED STATES	PART NUMBER	UPC NUMBER	CUFF SIZE
MASS MARKET	UA-787R	0-93764-60042-5	Medium 7.5" - 14.2" (19-36cm)
	UA-787V	0-93764-60026-5	Medium 7.5" - 14.2" (19-36cm)
HOME HEALTH CARE	UA-787H	0-93764-60027-2	Medium 7.5" - 14.2" (19-36cm)
CANADA	UA-787CN	0-93764-60028-9	Medium 7.5" - 14.2" (19-36cm)

Each unit comes with a trilingual instruction manual, quick start card, four week blood pressure logbook and warranty card. The UA-787R and UA-787V Mass Market models come with the same features and components, but are distinguished by different packaging. The UA-787H Home Health Care model also includes a carrying case and AC power adapter.

## Accessories

PART #	DESCRIPTION
UA-279	Small Blood Pressure Monitor Cuff: 5.1" - 7.9" (13-20cm)
UA-280	Medium Blood Pressure Monitor Cuff: 7.5" to 12.2" (19-31cm)
UA-281	Large Blood Pressure Monitor Cuff: 11.8" - 17.7" (30-45cm)
UA-380	<b>EasyCuff</b> feature: 7.5" - 14.2" (19-36cm)
UA-700CH	Cuff Hose
TB:181	AC Power Adapter

## Technical Specifications

Display	Digital: 20mm character height Pressure and pulse displayed simultaneously
Pressurization	Automatic, using micropump
Deflation	Active exhaust valve
Type	Oscillometric
Accuracy	Pressure: ±3mmHg or 2%, whichever is greater/Pulse: ±5%
Pressure Range	20mmHg to 280mmHg/Pulse: 40 to 200 pulses/minute
Power Source	Four type AA alkaline batteries (1.5 volt). Batteries not included.
Battery Life	Approx. four months with one daily measurement
Operating Environ.	50°F to 104°F (10°C to 40°C)
Storage Environ.	14°F to 140°F (-10°C to 60°C)

Blood pressure measurements determined with this device are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultation method within the limits prescribed by the American National Standards Institute (ANSI/AAMI SP10) for electronic sphygmomanometers. Manufactured by an ISO 9000 and ISO 13485 registered company. Engineered in Japan. Made in China.

## Packaging and Shipping Specifications

Selling Unit (Box)	L x W x H (inches)	WT (lbs.)	Individual Item Cube (cubic feet)
UA-787R/V/CN	5 <sup>1</sup> / <sub>8</sub> x 4 <sup>5</sup> / <sub>16</sub> x 7 <sup>5</sup> / <sub>8</sub>	1.40	0.10
UA-787H	6 <sup>15</sup> / <sub>16</sub> x 4 <sup>5</sup> / <sub>16</sub> x 7 <sup>11</sup> / <sub>16</sub>	2.20	0.13

Unit Case (10 units)	L x W x H (inches)	WT (lbs.)	Case Cube (cubic feet)
UA-787R/V/CN	22 x 11 x 8 <sup>1</sup> / <sub>2</sub>	15.50	1.19
UA-787H	26 x 15 x 11 <sup>1</sup> / <sub>4</sub>	25.90	2.54



1555 McCandless Drive  
Milpitas, CA 95035  
www.LifeSourceOnline.com

**LifeSource**  
Your source for a lifetime of health



LifeSource and EasyCuff are trademarks of A&D Medical. ©2002 A&D Medical. All rights reserved.



Ask your A&D Sales Representative about other **LifeSource** products.

- Blood Pressure Monitors
- Heart Rate Monitors
- Digital Thermometers
- Personal Scales
- Stethoscopes

For questions about **LifeSource** monitors, call our toll-free **LifeSource** HealthLine at 1-888-726-9966 or visit our web site at [www.LifeSourceOnline.com](http://www.LifeSourceOnline.com)

LIT:LS787 06/02